

Abutment Solutions

For customized implant restorations fabricated with CEREC® and inLab®

Digital all around.




**ivoclar
vivadent®**
passion vision innovation


The digital treatment workflow



Digital impression taking

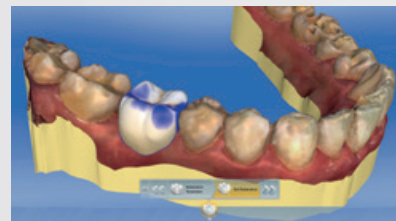

Digital planning

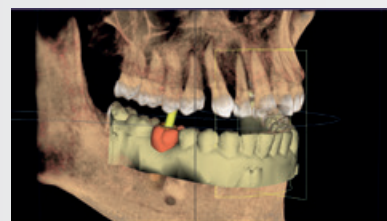

3D X-ray



Fabrication of the drilling template



Implantation*

Coordinated digital workflows in CAD/CAM technology open up new possibilities for implant-supported prosthetics – from impression taking to the final restoration – and ensure functional, highly esthetic, innovative and cost-efficient results. [1-3]



Implant design




Placement - Check-up*



Permanent restoration made of IPS e.max® CAD*



Shaped emergence profile following a temporary restoration*

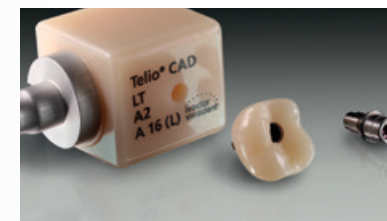

Temporary restoration made of Telio® CAD A16*

Together with Ivoclar Vivadent's innovative restorative materials, the CAD/CAM technique facilitates the planning and fabrication of implant-supported restorations – from temporary to permanent implant prosthetics. [2-4]

The following products are available for the fabrication of Abutment Solutions:

- Telio® CAD A16
- IPS e.max® CAD A14 / A16
- Multilink® Hybrid Abutment

*Source: Dr A. Kurbad [2], Germany

Dr S. Puri, USA



“ Telio CAD A16 and IPS e.max CAD blocks optimally complement each other. Therefore, CAD/CAM users are provided with a complete digital workflow which enables them to reliably fabricate temporary restorations and permanent hybrid abutment crowns. ”

Dr A. Kurbad, Germany



“ Telio CAD A16 blocks close the gap in the system chain of the proven IPS e.max CAD Abutment Solutions by including the temporization stage. The temporary restoration can be incorporated immediately after the implantation procedure or after the healing phase. Furthermore, it offers many options in terms of soft tissue management. Therefore, Telio CAD A16 forms the basis for an esthetic and functional therapy result. ”

The temporary –
key to a successful
treatment

Hybrid abutment crown made of Telio CAD A16:

- Cemented to the TiBase, suitable for immediate load-bearing or after the healing phase [2; 4]
- Easily designed emergence profile [2; 4]
- Visualization of the permanent restoration
- Blocks are available in the size A16 and in nine LT shades (BL3, A1, A2, A3, A3.5, B1, B3, C2, D2)



With Telio® CAD A16 from the temporary...

Telio® CAD A16 is designed for the CAD/CAM-supported fabrication of temporary hybrid abutment crowns. This 2-in-1 solution, i.e. crown and abutment combined, is used after the insertion of the implant and before the placement of the permanent restoration made of IPS e.max® CAD in the anterior and posterior region. [2; 4]

Telio CAD A16 supports the reconstruction of the gingiva during the first treatment phase. The surrounding soft tissue is individually formed. [2; 4]

The homogeneous, highly cross-linked PMMA block Telio CAD A16 features a predefined interface of size S or L for the direct cementation with the Dentsply Sirona TiBase.

Digital, individualized patient treatment in detail:




Clinical situation: Preparation for the digital, intraoral impression taking



Temporary hybrid abutment crown made of Telio® CAD A16




Seated Telio® CAD A16 restoration

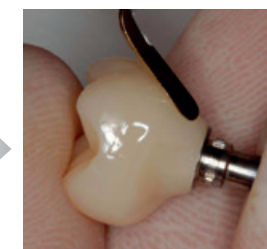


Optional: Design of the emergence profile by composite layering



Shaped emergence profile after removal of the temporary

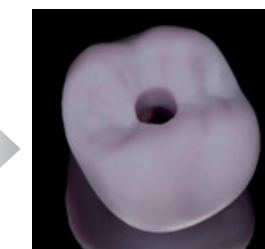

Milled abutment crown made of IPS e.max® CAD

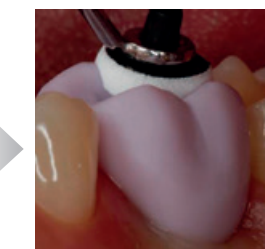

Clinical try-in



Final, seated IPS e.max® CAD hybrid abutment crown

Source: Dr L. Enggist / Dr Stephanie Huth, Ivoclar Vivadent AG

...to the final implant-supported IPS e.max® CAD restoration

Once a permanent solution is required, IPS e.max® CAD can be used to fabricate individual, implant-supported hybrid structures for single-tooth restorations using CAD/CAM technology. [1; 4; 5]

Depending on the indication, a tooth-shaded hybrid abutment with separate crown or a monolithic hybrid abutment crown can be fabricated. The particularly strong (530 MPa) [15] lithium disilicate glass-ceramic ties in seamlessly with the temporary Telio CAD restoration.

The IPS e.max CAD blocks A14 and A16 feature a predefined interface of size S or L for the direct cementation of the restoration with the Dentsply Sirona TiBase.


The permanent restoration –
flexibility during treatment

Hybrid abutment and hybrid abutment crown made of IPS e.max CAD:

- Excellent fit due to CAD/CAM processing technology [1-3]
- Esthetic restorations due to tooth-shaded hybrid abutment
- Hybrid abutment crown (2-in-1) offers functionality and efficiency [1; 4-8]
- Good biocompatibility with oral soft tissue [9-11]



Hybrid abutment

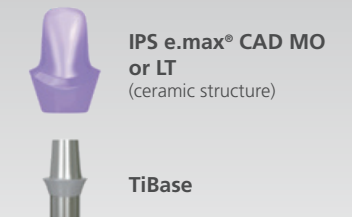

Hybrid abutment crown



Multilink® Hybrid Abutment – strong bond and esthetics



Telio® CAD and IPS e.max® CAD restorations are cemented to the TiBase in only a few steps.

Multilink® Hybrid Abutment

The self-curing luting composite Multilink Hybrid Abutment is used for the permanent cementation of ceramic and PMMA structures made of e.g. IPS e.max CAD or Telio CAD to TiBases.

This results in:

- a permanent cementation due to high bond strength values; [12; 13]
- optimal esthetics due to high opacity levels; [12]
- easy handling due to the convenient Automix syringe. [12]

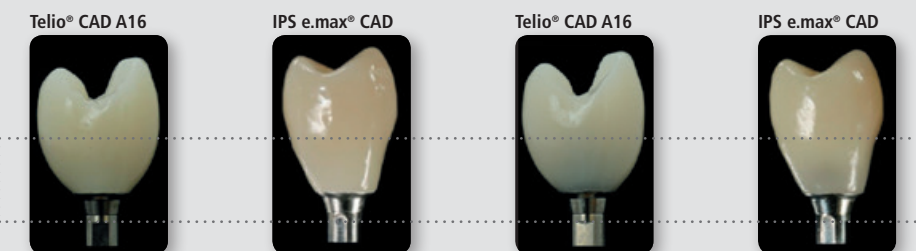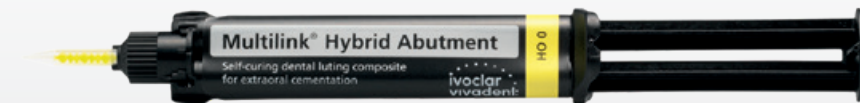

Optimum esthetics is achieved as the TiBase is entirely masked with **Multilink Hybrid Abutment**.

The opacity of other cementation materials is significantly too low.

SR Connect

The bonding agent conditions the temporary restoration made of Telio CAD A16 and prepares it for the cementation with Multilink Hybrid Abutment.

Monobond® Plus

The universal primer conditions the TiBase and IPS e.max CAD restoration and prepares both for the cementation with Multilink Hybrid Abutment.



IPS e.max® CAD Abutment Solutions – biocompatible and reliable


Physicochemical and biological study on the properties of polished versus glazed lithium disilicate ceramics (IPS e.max®)

University of Reims, France

The surfaces of the specimens made of IPS e.max (LS₂) glass-ceramic were given one of the following three surface finishes: untreated, polished or glazed. Thermanox served as control. The cell reaction on the polished and glazed surfaces was analyzed using a cell culture model based on chicken epithelium.

Result

Cell adhesion and proliferation (i.e. density) were higher on the polished than on the glazed surfaces. No cytotoxicity occurred in any of the specimens. Therefore, the high-strength lithium disilicate ceramic is not only a promising solution for esthetic implant abutments but it is also particularly suitable for sealing the periimplant bond.

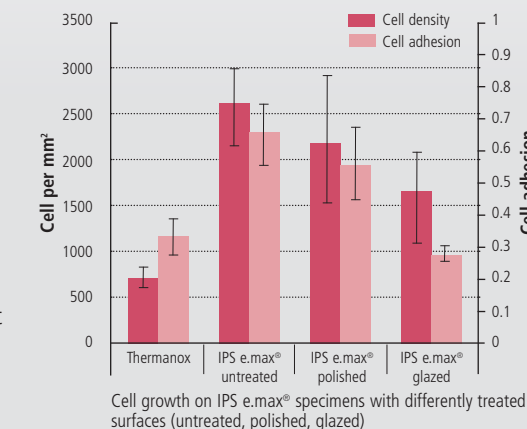

Source: C. Brunot-Gohin et al., see IPS e.max® Scientific Report Vol. 02/2001-2013

Study examining the reliability and failure types of ceramic abutments

New York University, USA

Specimens made of an IPS e.max (LS₂) hybrid abutment and screwed to an implant were prepared. IPS e.max (LS₂) crowns were cemented to the abutment.

Result

During the stress test in a universal testing machine, 100% of the hybrid abutments and hybrid abutment crowns made of IPS e.max (LS₂) glass-ceramic withstood a load of 280 N before fracture. The implant screw represented the system's weak point in all specimens. The screw fractured before any damage to the crown or abutment occurred.

Source: V.P. Thompson et al., [14], see IPS e.max® Scientific Report Vol. 02/2001-2013

References

1. Joda T, Bragger U. Time-efficiency analysis of the treatment with monolithic implant crowns in a digital workflow: a randomized controlled trial. Clin Oral Implants Res 2016;27:1401-1406.
2. Kurbad A. Provisional polymer-based CAD/CAM implant superstructures. CAD/CAM-basierte Implantatsuprastrukturen auf Polymerbasis als temporäre Versorgungen. Int J Comput Dent 2014;17:239-251.
3. Scherrer P, Watzke R, Enggist L, Heintze S, Tauch D. Designvalidierungstestbericht IPS e.max CAD A14/A16 Blocks for CEREC/inLab. Internal Investigation Report, Ivoclar Vivadent AG, Schaan, Principality of Liechtenstein, LL1753243 2013.
4. Çömlüoğlu ME, Nizam N, Çömlüoğlu MD. Immediate definitive individualized abutments reduce peri-implant bone loss: a randomized controlled split-mouth study on 16 patients. Clin Oral Investig 2017.
5. Joda T, Ferrari M, Bragger U. Monolithic implant-supported lithium disilicate (LS2) crowns in a complete digital workflow: A prospective clinical trial with a 2-year follow-up. Clin Implant Dent Relat Res 2017;19:505-511.
6. Zhang Y, Tian J, Wei D, Di P, Lin Y. Quantitative clinical adjustment analysis of posterior single implant crown in a chairside digital workflow: A randomized controlled trial. Clin Oral Implants Res 2019.
7. De Angelis P, Passarelli PC, Gasparini G, Boniello R, DiAmato G, De Angelis S. Monolithic CAD-CAM lithium disilicate versus monolithic CAD-CAM zirconia for single implant-supported posterior crowns using a digital workflow: A 3-year cross-sectional retrospective study. J Prosthet Dent 2019.
8. Mihali S, Canjau S, Bratu E, Wang H-L. Utilization of Ceramic Inlays for Sealing Implant Prostheses Screw Access Holes: A Case-Control Study. Int J Oral Maxillofac Implants 2016;31:1142-1149.
9. Fabel G, Beuer F. Establishing a suitable surface roughness for lithium disilicate implant abutments under laboratory conditions: a morphologic SEM and profilometric pilot study. Int J Comput Dent 2017;20:303-314.
10. Brunot-Gohin C, Duval JL, Azogui EE, Jannetta R, Pezron I, Laurent-Maquin D, Gangloff SC, Egles C. Soft tissue adhesion of polished versus glazed lithium disilicate ceramic for dental applications. Dent Mater 2013;29:e205-212.
11. Brunot-Gohin C, Duval J-L, Verbeke S, Belanger K, Pezron I, Kugel G, Laurent-Maquin D, Gangloff SC, Egles C. Biocompatibility study of lithium disilicate and zirconium oxide ceramics for esthetic dental abutments. Journal of periodontal & implant science 2016;46:362-371.
12. Watzke R, Heintze S, Albrecht T. Designvalidierungstestbericht Abutment Cement. Internal Investigation Report, Ivoclar-Vivadent AG, Schaan, Principality of Liechtenstein, LL1738501 2013.
13. Tauch D. Design Verifizierungsbericht – IPS e.max CAD A14 und A16. Internal Investigation Report, Ivoclar-Vivadent AG, Schaan, Principality of Liechtenstein, LL1773639 2013.
14. Thompson V, Coelho P, Silva N. Reliability and failure modes of a new ceramic abutment prototype. New York University, USA Research Report, Studienbericht für Ivoclar Vivadent, Data on File 2012.
15. Schweiger M. Report. Biaxial flexural strength of IPS e.max lithium disilicate products (CAD and Press). Internal Investigation Report, Ivoclar-Vivadent AG, Schaan, Principality of Liechtenstein, LL2943212 2016.

Ivoclar Vivadent AG
Bendererstr. 2
9494 Schaan
Liechtenstein
Tel. +423 235 35 35
Fax +423 235 33 60
www.ivoclarvivadent.com

673672/en/2019-11-28

