

Click & Bond

Adhese[®]
Universal

Top quality in
the efficient VivaPen[®]



Universal adhesion with advanced delivery

New user-friendly design

Adhese® Universal is a single-component, light-cured adhesive for direct and indirect bonding procedures and all etching protocols.^[1] It is the only universal bonding agent available in the efficient VivaPen® delivery form.

Fast and controlled intraoral application ^[2]

Click & bond:
Precise dispensing due to improved "click" mechanism

Unobstructed view of the operating field

Bendable tip gives easy access to all areas of the mouth and all cavity sizes

Exchangeable cannula with flocked tip

Fill level indicator

Your benefits at a glance:

- ✓ Universally suitable for direct and indirect bonding procedures and all etching protocols^[1]
- ✓ Consistent high bond strength (> 25 MPa) on wet and dry dentin
- ✓ Integrated desensitizing effect for minimizing^[3] the risk of postoperative sensitivity
- ✓ Efficient VivaPen: up to four times more applications per ml compared with conventional bottles


[1] Except in cases, where sufficient exposure to curing light cannot be ensured (e.g. adhesive cementation of endodontic posts).
[2] A protective sleeve must be used in intraoral procedures.
[3] Platt J.A. et. al., Dentaleconomics, 2016, p. 71-72

Quality that stands out

High bond strength for greater confidence

Adhese Universal achieves high bond strength values (> 25 MPa) on enamel and dentin irrespective of the etching protocol used or the moisture conditions on the tooth surface. Predictable, high bond strength values allow you to work with confidence.

Consistent high bond strength

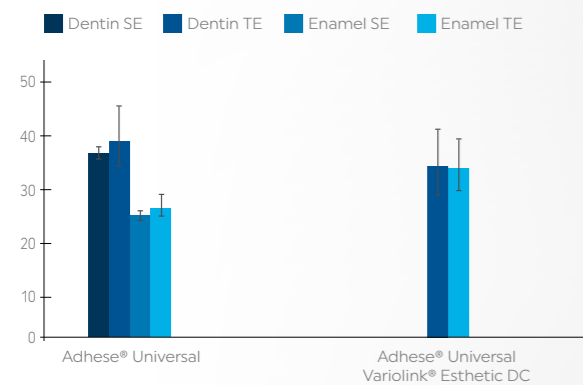

Direct bonding
Shear bond strength when using the self-etch (SE) technique and the total-etch (TE) technique
Source: DENTAL ADVISOR, USA, 2018

Indirect bonding – selective-enamel-etch technique
Shear bond strength with Variolink® Esthetic after dual-curing and 10,000 thermocycles (TC)
Source: R&D Ivoclar Vivadent AG, Schaan, 2015


Minimizing the risks of postoperative sensitivity^[3]

Adhese Universal features an integrated desensitizing effect, which eliminates the need for a separate desensitizing agent. Adhese Universal forms a mechanical barrier and seals the dentin tubules under wet and dry conditions.

Very low film thickness for fitting accuracy and esthetic results

Since Adhese Universal can be applied very thinly, you have all the freedom you need to create esthetic direct and indirect restorations. The adhesive establishes a sound foundation for achieving customized solutions and esthetic results.

Shear bond strength on dentin after different surface treatments – total-etch technique¹

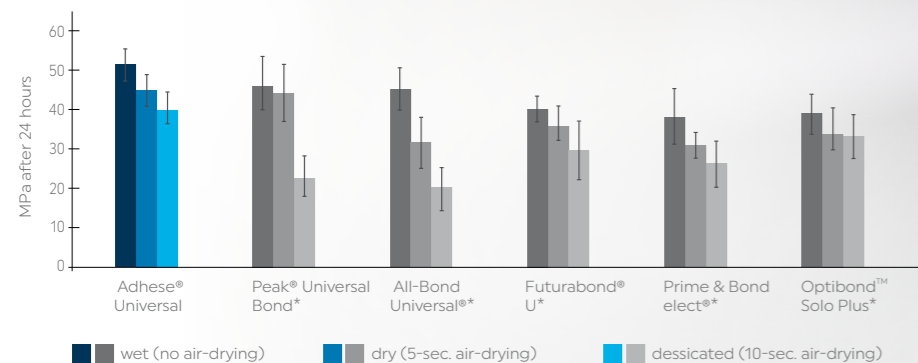

* Not registered trademarks of Ivoclar Vivadent AG. Source: Singhal S et al., Dent. Mater. 2014, 30, p. e81-e82.



Remember:

Before light curing, dry with oil- and water-free, clean compressed air until the etched fusion surfaces are chalk-white!


[3] Platt J.A. et. al., Dentaleconomics, 2016, p. 71-72

Threefold efficiency

Up to
4x more
applications
[4-6]

Economically dispensed with VivaPen

The simple "click" mechanism pumps the exact amount of adhesive needed. The adhesive does not have to be predispensed into a mixing well. This significantly reduces the amount of material wasted. As a result, up to four times more single tooth applications are obtained per ml of pen contents compared to conventional bottles.

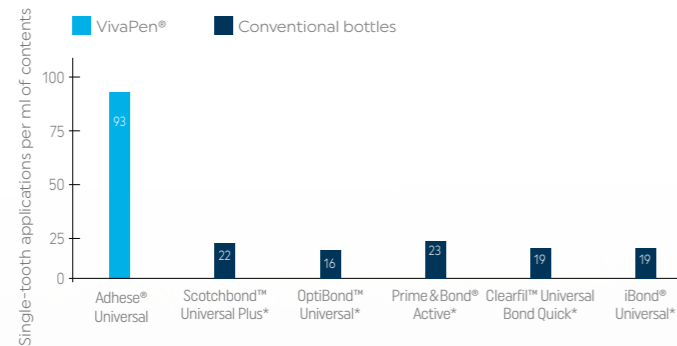

The measurements were taken using gravimetric methods and simulating daily clinical use.
* Not registered trademarks of Ivoclar Vivadent AG.
Sources: Berndt & Partner, VivaPen Benchmarking Study, Study Report, April 2022.

Universal
use

For direct and indirect restorations and all etching protocols^[7]

Adhese Universal can be used in conjunction with light-curing, self-curing and dual-curing composite restorative materials. It is the adhesive of choice for bonding composites of the Tetric line. The product acts as an adhesive in the placement of indirect and direct restorations with light and dual-curing luting materials or as a primer for conditioning restorations made from Tetric CAD composite blocks.

3-second
cure

Fast, reliable cure

Adhese Universal is cured with conventional curing lights (with a light intensity of 2000 mW/cm²) in just 5 seconds. For restorations of Class I and II cavities in permanent posterior teeth with occlusal ligh curing using the Bluephase PowerCure with a light intensity of 3000 mW/cm² the curing process is reduced to only 3 seconds.

Learn more about our 3s PowerCure system of products:
ivoclar.com/3sPowerCure

Adhese Universal DC

Strong in any situation

Adhese Universal DC is a single-component, dual-cured adhesive for universal use. It is particularly suitable for bonding indirect restorations in situations where self-curing is required or desired.^[8]

Same bonding performance as Adhese Universal

- ✓ Available in Free Stand® Single Dose units with initiator-coated applicators in two sizes.
- ✓ Even better together – coordinated with Variolink Esthetic for the placing of all types of restorations and root canal posts. Can be used in the self-cure mode with Variolink Esthetic DC and all etching protocols.



Delivery forms

Adhese Universal	
754798WW	Kit VivaPen (1 x 2 ml) 100 VivaPen snap-on cannulas, 20 VivaPen protective sleeves
754799WW	Starter Kit VivaPen (1 x 0.5 ml) 10 VivaPen snap-on cannulas, 20 VivaPen protective sleeves
Adhese Universal DC	
687778WW	Adhese Universal DC Refill Single Dose Regular (25 x 0.1g)
687779WW	Adhese Universal DC Refill Single Dose Small/Endo (25 x 0.1g)

Additional delivery forms are listed here: ivoclar.com/AdheseUniversal/deliveryforms

[4] per ml compared with conventional bottles

[5] Berndt & Partner, Benchmarking Study, Study Report, April 2022.

[6] Holaschke S, VEA Verification Report, Ivoclar Vivadent, 2021.

[7] Except in cases, where sufficient exposure to curing light cannot be ensured (e.g. adhesive cementation of endodontic posts).

[8] Compatible with light and dual-curing luting composites using an etch-and-rinse protocol.

Direct restorative procedures for efficient esthetics

Harness the power of our coordinated products for direct restorative procedures to create high-quality esthetic restorations efficiently.



Isolate

Comfortable relative isolation with OptraGate



Bond

Efficient bonding with the Adhese Universal VivaPen



Restore

Tetric Line – The solution for all cavity classes



Contour

Effortless contouring with OptraSculpt



Cure

Powerful performance with Bluephase PowerCure



Polish

Shiny surfaces with OptraGloss



Protect

Effective prevention with Fluor Protector S

The range of available products may vary from country to country.

Ivoclar Vivadent AG | Bendererstr. 2 | 9494 Schaan | Liechtenstein
Tel. +423 235 35 35 | Fax +423 235 33 60

ivoclar.com

666987EN/2023-10-03