

**INTRODUCING**

# Tetric<sup>®</sup> CAD

Esthetic Hybrid Ceramic Resin

**Inlay & Onlay,  
*the faster way!***

**DIGITAL EXPERTISE  
BEGINS WITH MATERIALS**



**INTRODUCING Tetric<sup>®</sup> CAD Hybrid Ceramic Resin block**

- **EFFICIENT** – fast milling and polishing
- **ESTHETIC** – medium & high translucency for superior esthetics
- **HIGH STRENGTH** – flexural strength 272 MPa
- **PROVEN CEMENTATION** – with Adhese<sup>®</sup> Universal and Variolink<sup>®</sup> Esthetic

**ivoclar**  
**vivadent**<sup>®</sup>  
passion vision innovation


# Tetric<sup>®</sup> CAD

Esthetic Hybrid Ceramic Resin Block

Tetric<sup>®</sup> CAD is a hybrid ceramic resin block for the use of efficient fabrication of single-tooth restorations. Due to the pronounced chameleon effect, Tetric CAD restorations blend in beautifully with the surrounding tooth structure. Restorations are milled, polished and then seated using an adhesive cementation technique. This processing procedure is very efficient, leading to quick and easy esthetic results.


## Flexural Strength

Comparison of biaxial flexural strength of Tetric CAD and various hybrid materials to the standard ISO 6872. Tetric CAD exhibited a high flexural strength of 272 MPa.

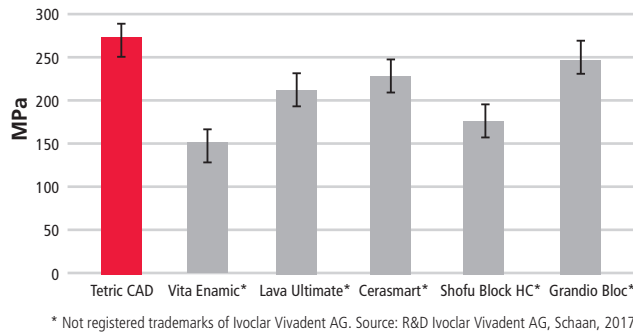

## Water Absorption

The water absorption of Tetric CAD and various hybrid materials were tested according to the standard ISO 10477. Water absorption was well below the level specified by the standard (<40µg / mm<sup>3</sup>), for most of the hybrid materials including Tetric CAD. Lava Ultimate exhibited the highest water absorption levels in this study.

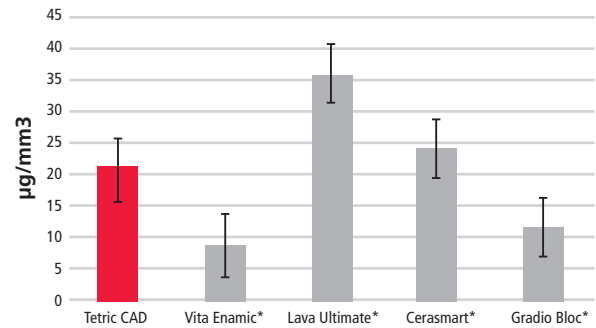

**Indications** Veneers, inlays, onlays (e.g. partial crowns, occlusal veneers), crowns

### HT (High Translucency)




### Shades



### MT (Medium Translucency)



### Shades

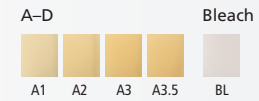

## Tetric<sup>®</sup> CAD | Introductory Specials

**BUY** 4 Refills of Tetric<sup>®</sup> CAD for CEREC<sup>®</sup>

**BUY** 4 Refills of IPS e.max<sup>®</sup> CAD (LT, HT, MT)

**BUY** 4 Refills of Tetric<sup>®</sup> CAD for CEREC<sup>®</sup>

**BUY** 4 Refills of Tetric<sup>®</sup> EvoCeram<sup>\*</sup> or Tetric<sup>®</sup> EvoCeram Bulk Fill<sup>\*</sup>

\*Refills must be same brand and delivery form. Can mix and match shades.

**BUY** 1 Variolink<sup>®</sup> Esthetic Promo Pack

**BUY** 3 Tetric<sup>®</sup> CAD Refills for CEREC<sup>®</sup>

**GET** 1 Refill of Tetric<sup>®</sup> CAD for CEREC<sup>®</sup>

**GET** 1 Refill of Tetric<sup>®</sup> CAD for CEREC<sup>®</sup>

**GET** 1 Adhese<sup>®</sup> Universal System (VivaPen<sup>®</sup>)

**GET** 1 Refill of Tetric<sup>®</sup> CAD for CEREC<sup>®</sup>

**GET** 1 Refill of Tetric<sup>®</sup> CAD for CEREC<sup>®</sup>

**GET** 1 Variolink<sup>®</sup> Esthetic 5g Refill

### FREE GOOD REDEMPTION: (VALID IN USA ONLY)

All purchases must appear on a single invoice between 8/1/18-12/30/18. Requests for free goods must be received by 1/31/19, by forwarding a copy of your qualifying invoice to: Ivoclar Vivadent, Inc., Attn: 2018 Clinical Product Specials, 175 Pineview Drive, Amherst, NY 14228, FAX: (800) 598-4569 e.mail: usclinical@ivoclarvivadent.com. Ivoclar Vivadent reserves the right to modify or cancel these promotions at any time. These specials cannot be combined with any other offer or rebate program.

For more information, call us at 1-800-533-6825 in the U.S., 1-800-263-8182 in Canada.

©2018 Ivoclar Vivadent, Inc. Variolink, Tetric, Adhese, VivaPen, IPS e.max, are all trademarks of Ivoclar Vivadent, Inc. CEREC is not a registered trademark of Ivoclar Vivadent, Inc. 701257 Rev. 1 8/2018

